


Model: CMRTSS400

## 400L Solid Door RTS Cabinet

### Overview

- +15°C to +25°C temperature range
- Core temperature display, represents your vaccine vial through **iCOOL** INTELLIGENT COOLING temperature probe
- External digital controller and temperature display
- Battery back up
- Annual service & Calibration available
- Fitted lock with 2 keys

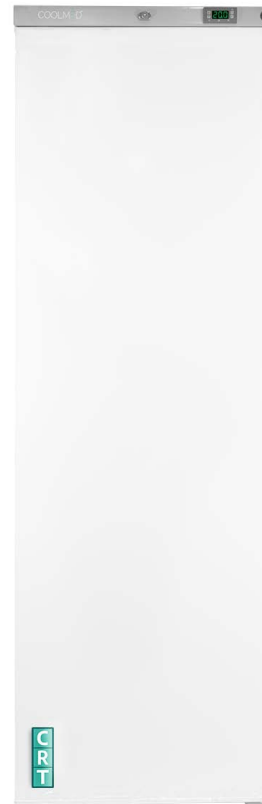

**2 YEAR PARTS & LABOUR WARRANTY**

**ANNUAL SERVICE & CALIBRATION AVAILABLE**

### Key Features




#### Controller

Controller system. Read, record, reset with the touch of one button. Adjustable temperature setting.



#### Battery back up

Battery back up should there ever be a power outage the battery backup will continue to monitor the internal temperature of the fridge.

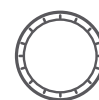

#### Porthole

Factory fitted porthole enables quick and easy access for any secondary temperature monitoring probes.



#### LED lighting

Low voltage LED lighting improves running cost as no heat is produced. Long lasting.



#### Flexible storage

Multiple shelf locations. Back stop on every shelf provides a flexible storage solution.

## Details

Model: CMRTSS400

Gross Capacity (litres)	408	Number of Shelves	7 (6 full size, 1 small base shelf)
Net Capacity (litres)	400	Number of Shelf Positions	10+
Unit External Width (mm)	595	Shelf Material	Plastic Coated Wire
Unit External Depth (mm)	695	Interior Light	Yes
Unit External Height (mm)	1860	Light Type	LED
Product Weight (kg)	74	Light Activation	Switch
Packed Width (mm)	650	Fan	Yes
Packed Depth (mm)	750	Spark Free	No
Packed Height (mm)	1921	Controller Type	Digital, External
Packed Weight (kg)	82	Number of Probes	2
Temperature Range	+15°C to +25°C	Probe Type	Air & Load
External Finish, Internal Finish	Painted Steel, Painted Steel	Ambient Temperature Range	+10°C to +30°C
Battery Back Up	Yes	Refrigerment Type	R600a
Door Finish, Door Type, Door Style	Painted Steel, Solid, Hinged, Right Hand	CFC Free, HCFC Free, Ammonia Free	Yes
Open Door Sensor	Yes	Voltage Supply	220-240V/ 50Hz
Lock	Yes	Certification	CE/UK
Adjustable Feet (number)	Yes (2)	Noise Level (dBA)	50
Porthole	Yes	Input Power (W)	240
Porthole Location	Rear, Middle	Cable Length (cm)	20
Temperature Display on Panel, Temperature Alarm, Remote Alarm Contacts, Audible Alarms, Visual Alarms	Yes	Energy Consumption (kWh/annum)	2.6

## Related Products



**CMRTSG400**  
400L Glass Door  
RTS Cabinet



**CMDL01**  
Advanced Data Logger  
with LCD Display



**COMBI01**  
Keyless  
Combination Lock