

Model: CMRTSG300

## 300L Glass Door RTS Cabinet

### Overview

- +15°C to +25°C temperature range
- Core temperature display, represents your vaccine vial through **iCOOL** INTELLIGENT COOLING temperature probe
- External digital controller and temperature display
- Battery back up
- Annual service & Calibration available
- Fitted lock with 2 keys

**2 YEAR PARTS & LABOUR WARRANTY**

**ANNUAL SERVICE & CALIBRATION AVAILABLE**



### Key Features




#### Controller

Controller system. Read, record, reset with the touch of one button. Adjustable temperature setting.



#### Battery back up

Battery back up should there ever be a power outage the battery backup will continue to monitor the internal temperature of the fridge.

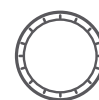

#### Porthole

Factory fitted porthole enables quick and easy access for any secondary temperature monitoring probes.



#### LED lighting

Low voltage LED lighting improves running cost as no heat is produced. Long lasting.



#### Flexible storage

Multiple shelf locations. Back stop on every shelf provides a flexible storage solution.

## Details

Model: CMRTSG300

Gross Capacity (litres)	308
Net Capacity (litres)	300
Unit External Width (mm)	595
Unit External Depth (mm)	650
Unit External Height (mm)	1495
Product Weight (kg)	72
Packed Width (mm)	650
Packed Depth (mm)	702
Packed Height (mm)	1556
Packed Weight (kg)	78
Temperature Range	+15°C to +25°C
External Finish, Internal Finish	Painted Steel, Painted Steel
Battery Back Up	Yes
Door Finish, Door Type, Door Style	Painted Steel, Glass, Hinged, Right Hand
Open Door Sensor	Yes
Lock	Yes
Adjustable Feet (number)	Yes (2)
Porthole	Yes
Porthole Location	Rear, Middle
Temperature Display on Panel, Temperature Alarm, Remote Alarm Contacts, Audible Alarms, Visual Alarms	Yes

Number of Shelves	7 (6 full size, 1 small base shelf)
Number of Shelf Positions	10+
Shelf Material	Plastic Coated Wire
Interior Light	Yes
Light Type	LED
Light Activation	Switch
Fan	Yes
Spark Free	No
Controller Type	Digital, External
Number of Probes	2
Probe Type	Air & Load
Ambient Temperature Range	+10°C to +30°C
Refrigerment Type	R600a
CFC Free, HCFC Free, Ammonia Free	Yes
Voltage Supply	220-240V/ 50Hz
Certification	CE/UK
Noise Level (dBA)	50
Input Power (W)	220
Cable Length (cm)	20
Energy Consumption (kWh/annum)	2.2

## Related Products



**CMRTSS300**  
300L Solid Door  
RTS Cabinet



**CMDL01**  
Advanced Data Logger  
with LCD Display



**COMBI01**  
Keyless  
Combination Lock