




Model: CMS29

29L Small Solid Door Refrigerator

Overview

- +2°C to +8°C temperature range
- Core temperature display, represents your vaccine vial through **iCOOL** INTELLIGENT COOLING temperature probe
- External digital controller and temperature display
- Battery back up
- Self-closing door
- Annual service & calibration available
- Fitted lock with 2 keys

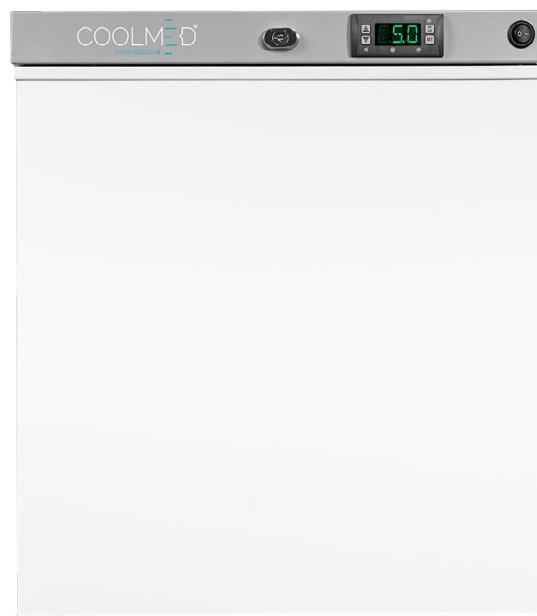

5 YEARS PARTS & LABOUR WARRANTY

ANNUAL SERVICE & CALIBRATION AVAILABLE

Key Features



Controller

Easy-to-use control system with a one-touch button to read, record, and reset settings. Includes fully adjustable temperature controls for precise monitoring.



Battery back up

In the event of a power outage, the integrated battery backup ensures continuous monitoring of the internal temperature—keeping your contents safe and secure.



Porthole

Pre-fitted factory porthole allows simple, hassle-free access for adding secondary temperature monitoring probes.



LED lighting

Energy-efficient, low-voltage LED lighting reduces running costs and eliminates excess heat. Durable and long-lasting for everyday use.



Flexible storage

Customise your space with multiple shelf positions. Each shelf includes a back stop to prevent items from being pushed too far—practical and versatile.

Details

Model: CMS29

Gross Capacity (litres)	30
Net Capacity (litres)	29
Unit External Width (mm)	448
Unit External Depth (mm)	484
Unit External Height (mm)	535
Product Weight (kg)	23
Packed Width (mm)	490
Packed Depth (mm)	510
Packed Height (mm)	580
Packed Weight (kg)	25
Temperature Range	+2°C to +8°C
External Finish, Internal Finish	Painted Steel, Painted Steel
Battery Back Up	Yes
Door Finish, Door Type, Door Style	Painted Steel, Solid, Hinged, Right Hand
Open Door Sensor	Yes
Lock	Yes
Adjustable Feet (number)	Yes (2)
Porthole	Yes
Porthole Location	Rear, Middle
Temperature Display on Panel, Temperature Alarm, Remote Alarm Contacts, Audible Alarms, Visual Alarms	Yes

Number of Shelves	3 (2 full size, 1 small)
Number of Shelf Positions	6
Shelf Material	Plastic Coated Wire
Interior Light	Yes
Light Type	LED
Light Activation	Switch
Fan	Yes
Spark Free	No
Controller Type	Digital, External
Number of Probes	2
Probe Type	Air & Load
Ambient Temperature Range	+10°C to +30°C
Refrigerant Type	R600a
CFC Free, HCFC Free, Ammonia Free	Yes
Voltage Supply	220-240V/ 50Hz
Certification	CE/UK
Noise Level (dBA)	50
Input Power (W)	76
Cable Length (cm)	20
Energy Consumption (kWh/annum)	0.6

Related Products

**CMG29**Small Glass
Door Refrigerator**CMDL01**Advanced Data Logger
with LCD Display**COMBI01**Keyless
Combination Lock
www.coolmed.co.uk
+44 161 772 5666 | sales@coolmed.co.uk

8 Coronet Way, Salford, Manchester, M50 1RE

coolmed.co.uk